



## Autoclaves for sterilization "MICRO 8", "MED 12" and "MED 20"

PRESSURE REGULATOR THERMOSTAT OF 1 AND 2 bar. TEMPERATURES 121 °C AND 134 °C.  
 VERTICAL DESKTOP MODELS. "MICRO 8" CAPACITY: 8 LITRES. "MED 12" CAPACITY: 12 LITRES.  
 VERTICAL STAND MODEL "MED 20" CAPACITY: 20 LITRES.

*An economic alternative. Semiautomatics.*

### SAFETY:

**SAFETY VALVE:** PREVENTS THE PRESSURE WITHIN THE AUTOCLAVE OVERPASSING THE MAXIMUM SET LIMIT.  
**SAFETY PRESSURE:** DISCONNECTS POWER TO THE HEATER IN CASE OF OVERPASSING THE PRESSURE SET LIMIT.  
**SAFETY THERMOSTAT:** DISCONNECTS THE HEATING IN CASE OF OVERTEMPERATURE OR LACK OF WATER IN THE TANK.  
 THERMALLY INSULATED LID MADE OF MOLDED THERMOPLASTIC.

### COMMON FEATURES

AISI 304 stainless steel external case and top.  
 Stainless steel chamber and lid, with easy opening.  
 INCOLOY special stainless steel heater element, which provides long life, low corrosion and high temperatures resistance.  
 Manual air purge when starting ("MICRO 8" model).  
 Automatic atmospheric purge ("MED 12" and "MED 20" models).  
 Safety device to prevent lid opening.  
 Starting main switch.




#### REAR PART

Mains source.  
 Two fuse carriers.  
 Safety thermostat.

#### SIDE PART

Drain outlet.  
 Steam and drainage selective valve.

### CONTROL PANEL

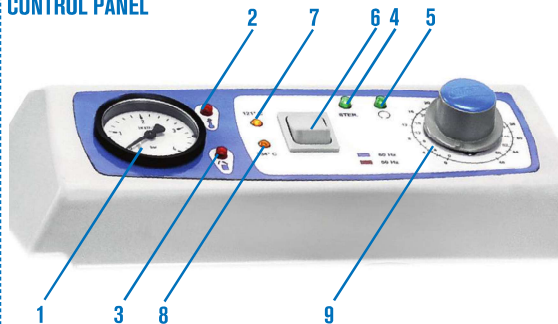

1. Pressure gauge.
2. Over temperature alarm light indicator.
3. Door open light indicator.
4. Sterilization phase light indicator.
5. Sterilization end of cycle light indicator.
6. Temperature selector 121°C or 134°C.
7. Light indicator at 121°C.
8. Light indicator at 134°C.
9. Timing from 0 to 60'.



## Autoclave for sterilization "MICRO 8"

DESKTOP MODEL



### MODEL

Part No.	Capacity litres	Ø (usable) cm	Height (usable) cm	Height/Width/Depth (exterior) cm			Power W	Weight Kg
4001745	8	24	15	40	38	38	1000	13

It comes complete with a stainless steel wire basket of Ø 22x14cm high. Part no. 1004771.



## Autoclave for sterilization "MED 12"

DESKTOP MODEL

Features: see page 56



### MODEL

Part No.	Capacity litres	Ø (usable) cm	Height (usable) cm	Height/Width/Depth (exterior) cm	Power W	Weight Kg
<b>4001756</b>	12	24	26	46 38 38	1800	18

It comes complete with a stainless steel wire basket of Ø 22x24cm high. Part no. 1004772.




## Autoclave for sterilization "MED 20"

VERTICAL STAND MODEL

CONFORMS TO APPARATUS PRESSURE REGULATION DIRECTIVES.

Features: see page 56

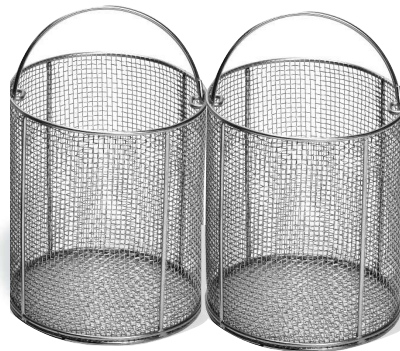

### MODEL

Part No.	Capacity litres	Ø (usable) cm	Height (usable) cm	Height/Width/Depth (exterior) cm	Power W	Weight Kg
<b>4001757</b>	20	24	41	70 38 38	2500	41

It comes complete with two stainless steel wire baskets of Ø 23x17 cm Part No. 1004773