


Advantages of the Smart controller:

- 4,3", clear, full colour touch screen
- USB and LAN ports for data transfer
- multi-segment time and temperature programs
- loop function up to 255 times or endless
- adjustable start delay feature
- adjustable ramps
- recording of min, average and max temperature value for each segment
- overview of set and current parameters while operating
- operating in temperature or time priority mode
- fan speed control in the range of 0 ... 100% or 10% ... 100% (depending on the chamber capacity) except ST/CHL/ZL
- temperature sensor fail alarm
- possibility of temperature calibration by the User
- power failure control system (program continued after restoring power)
- visual and sound alarm
- internal memory for programs and data storage
- operating with gloves on
- event registry
- user manual for direct download
- ability to buy a license for [LabDesk](#) software
- Quick change of program parameters
- Alarm Bar



Touch screens of the Smart and Smart PRO controllers can be operated with latex gloves!

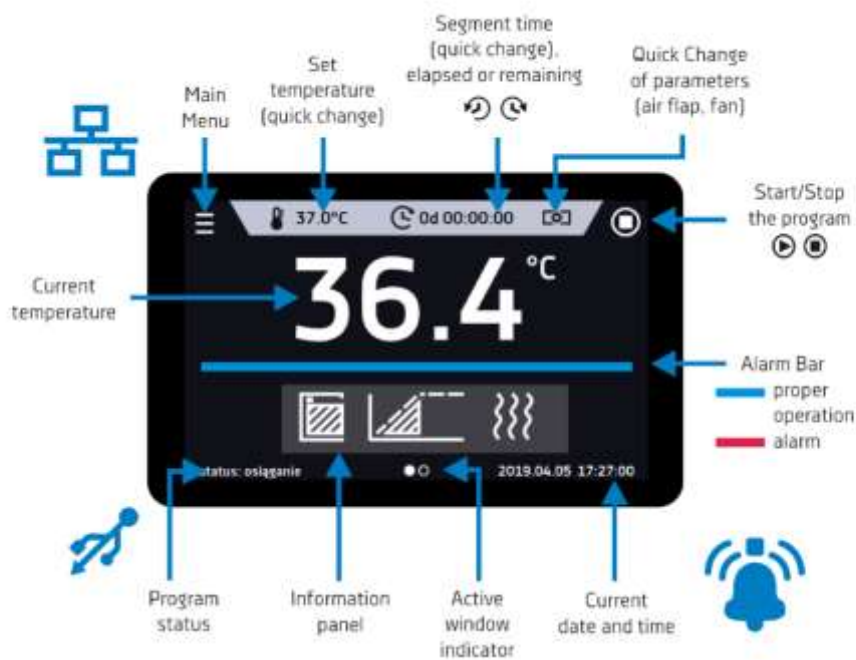

The Smart and Smart PRO controllers feature wide angle displays. You can see the screen clearly from any angle, regardless of the light intensity.