

## Manual pipette controller


### The Profiller™ manual 437

**Robust and easy to use, this model offers ergonomic shape and lightweight. It provides smooth control in both aspirating and dispensing, using glass or plastic pipettes of volume up to 100 mL. One-year warranty.**

- Large capacity bulb
- Soft action lever button
- Efficient blow-out system
- PTFE protection membrane filter
- Maintenance free

**437**

Pipettes  
up to 100 mL

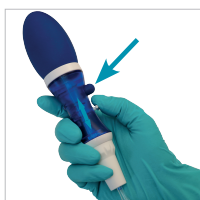

#### Squeezable bulb

Simply squeeze large capacity silicone bulb for efficient aspiration. If additional aspiration is needed, it can be squeezed even during pipetting.

#### Precise control

Thumb lever button for left and right-handed users to easily control both aspiration and dispensing speeds.

#### Liquid blow-out

Press large button to expel residual liquid when using blow-out type pipettes.

#### Protection filter

Interchangeable hydrophobic membrane filter guarantees effective protection against airborne contamination during aspiration and dispensing.

#### Pipette holder

Conical shape of silicone pipette holder tightly fits graduated glass or plastic pipettes. Autoclavable.

### Ordering information - Profiller™ manual 437

| Description                      | Packaging | Cat. No. |
|----------------------------------|-----------|----------|
| Profiller™ manual 437, blue bulb | 1 / pk    | 437100   |

### The Profillette™ 406

**Low cost, easy to use pipette filler suitable for graduated glass and plastic pipettes. Rotate wheel to aspirate. Fast release by pressing air inlet valve button, or gradual release by gently turning release wheel.**



**406**

2 mL  
10 mL  
25 mL

- One-hand operation
- Easy filling and dispensing
- Air inlet valve for fast distribution
- Colour coding per size
- Maintenance free

### Ordering information - Profillette™ 406

| Volume         | Colour   | Packaging | Cat. No. |
|----------------|----------|-----------|----------|
| 2 mL           | blue     | 1 / pk    | 406.002  |
| 10 mL          | green    | 1 / pk    | 406.010  |
| 25 mL          | red      | 1 / pk    | 406.025  |
| 2 / 10 / 25 mL | assorted | 3 / pk    | 406.300  |

### Ordering information - Accessories

| Description                   | Packaging | Cat. No. |
|-------------------------------|-----------|----------|
| <b>Profillette™ 406</b>       |           |          |
| Spare silicone pipette holder | 1 / pk    | 1.406.01 |
| Replacement nose screw        | 1 / pk    | 1.406.02 |
| <b>Profiller™ manual 437</b>  |           |          |
| Membrane filter, 0.45 µm      | 5 / pk    | 322.447  |
| Silicone pipette holder       | 1 / pk    | 1.437.01 |
| Nose cone - white             | 1 / pk    | 1.437.02 |