

Sterilization pouches


The Qualitix® sterilization pouches

Fitting any small to mid-size laboratory, medical, dental or veterinary device, the Qualitix® single use pouches provide easy and effective wrapping of items requiring sterilization. Made with robust medical grade paper, the pouches tightly close with self-sealing strip to maintain sterility. Content remains always visible through PET / CPP transparent laminated film.

- Robust durable material
- Quick content identification
- Large self-sealing strip, easy to close
- Chevron top seal for easy opening
- Suitable for steam and EO gas sterilization
- Two sterilization indicators
- Five different sizes
- Conformity to ISO 17665, ISO 11135, MDR 2017/745



Free samples on request

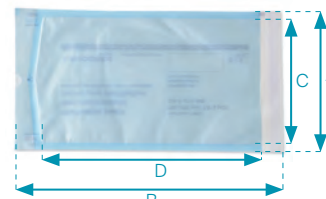

Sterilization indicators

Two chemical indicators included to confirm proper sterilization process. Blue indicator turns dark grey after autoclaving, red indicator turns to yellow after ethylene oxide (EO) gas sterilization.

Integrity control

Blue stained laminated film enables easy control of film/pouch integrity. Do not use a damaged pouch.

Self-sealing - quick opening

Simply peel off protection strip, fold and press with fingers. Pouch remains tightly sealed to maintain product sterility.

Dimensions

A (mm)	B (mm)	C (mm)	D (mm)
83	160	65	110
90	260	70	205
135	260	115	205
190	360	170	300
250	370	220	275
300	395	270	335
300	490	270	390

Ordering information

Description	Dimensions		Material	Packaging	Cat. No.
Pouch, self-sealing	83 x 160 mm	3.3" x 6.3"	Med. grade paper/PET/CPP	200 / pk	3210.0816
Pouch, self-sealing	90 x 260 mm	3.5" x 10.2"	Med. grade paper/PET/CPP	200 / pk	3210.0926
Pouch, self-sealing	135 x 260 mm	5.3" x 10.2"	Med. grade paper/PET/CPP	200 / pk	3210.1326
Pouch, self-sealing	190 x 360 mm	7.5" x 14.2"	Med. grade paper/PET/CPP	200 / pk	3210.1936
Pouch, self-sealing	250 x 370 mm	9.8" x 14.6"	Med. grade paper/PET/CPP	200 / pk	3210.2537
Pouch, self-sealing	300 x 395 mm	11.8" x 15.6"	Med. grade paper/PET/CPP	200 / pk	3210.3040
Pouch, self-sealing	300 x 490 mm	11.8" x 19.3"	Med. grade paper/PET/CPP	200 / pk	3210.3049



See video